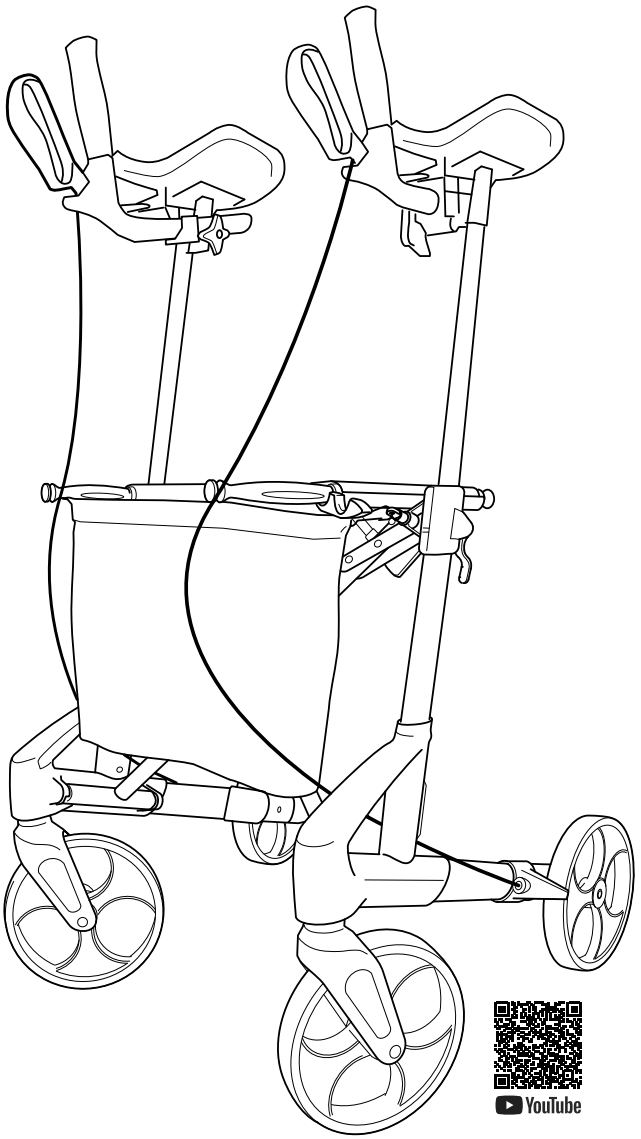
EN

SV

TOPRO

# Instructions for use
# TOPRO Troja Walker² (SE)

7 YEAR WARRANTY

CE

topromobility.com

made in Norway

EN

# Illustrations
# Dimensions
# Technical data
# Article numbers

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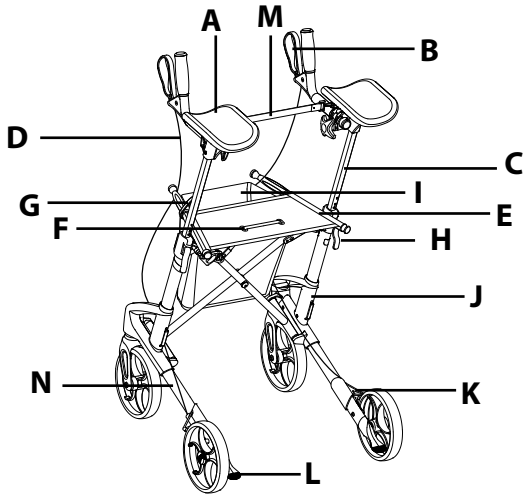
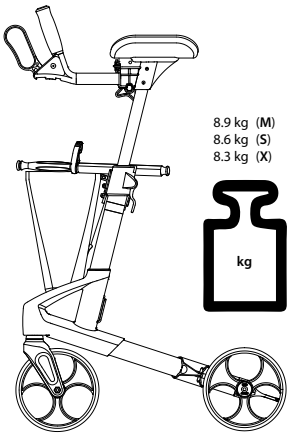
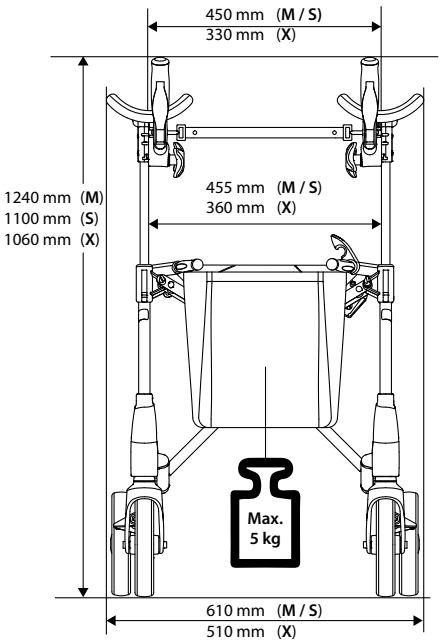
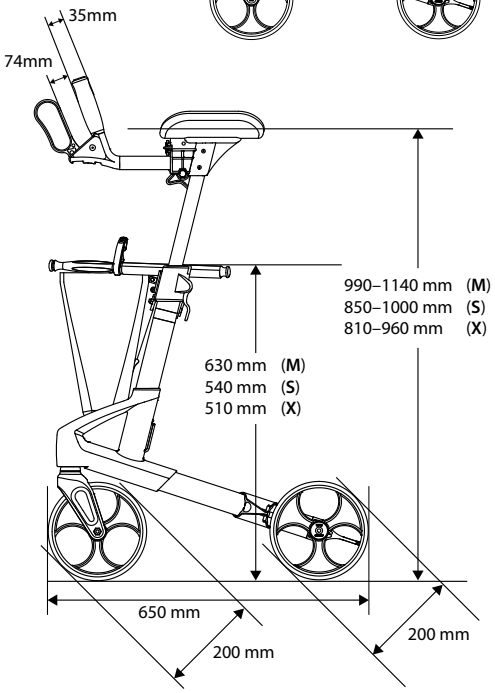
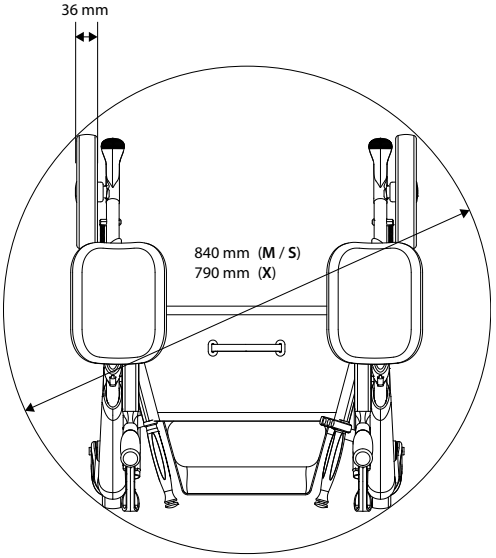
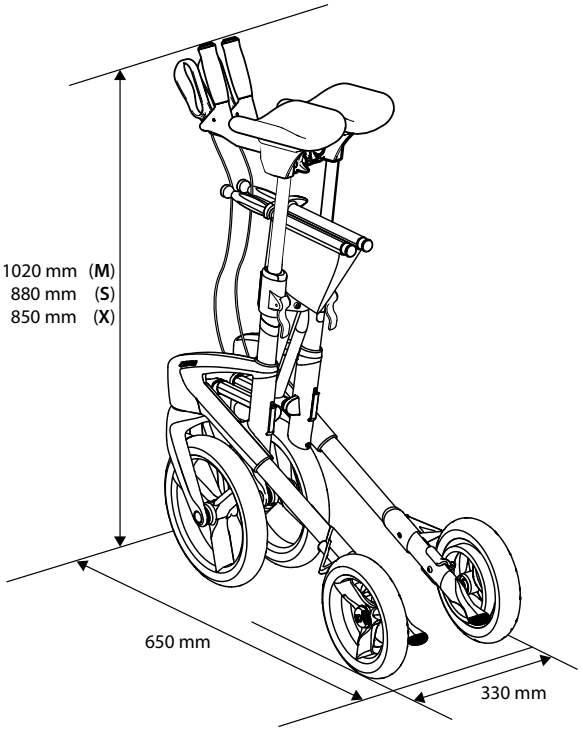
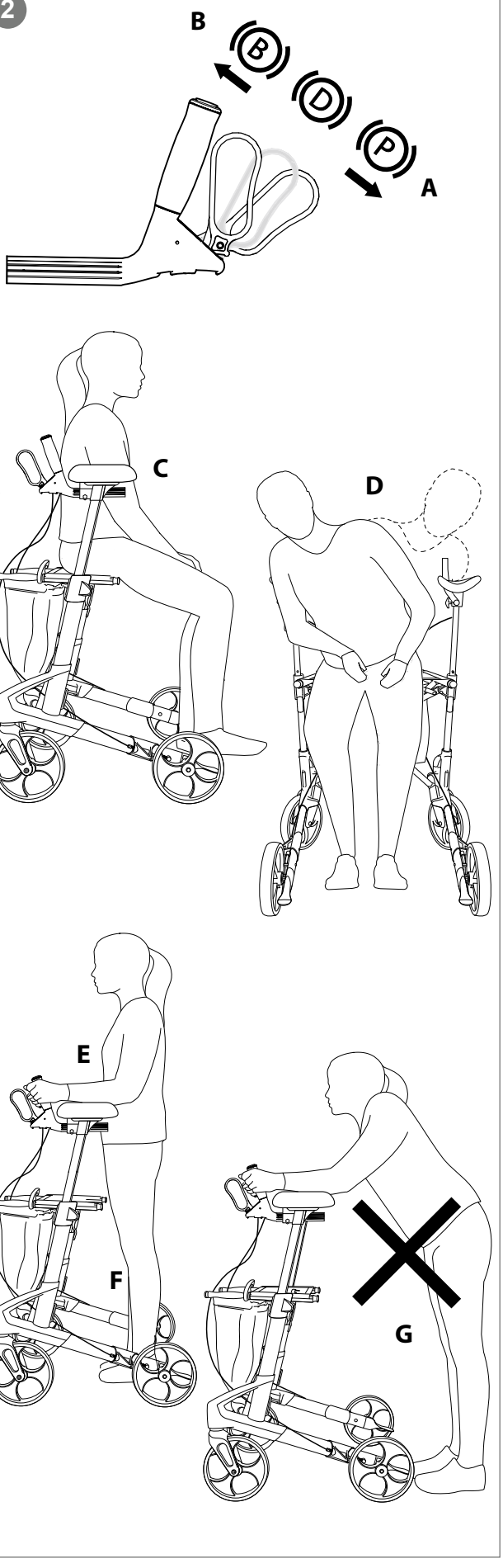
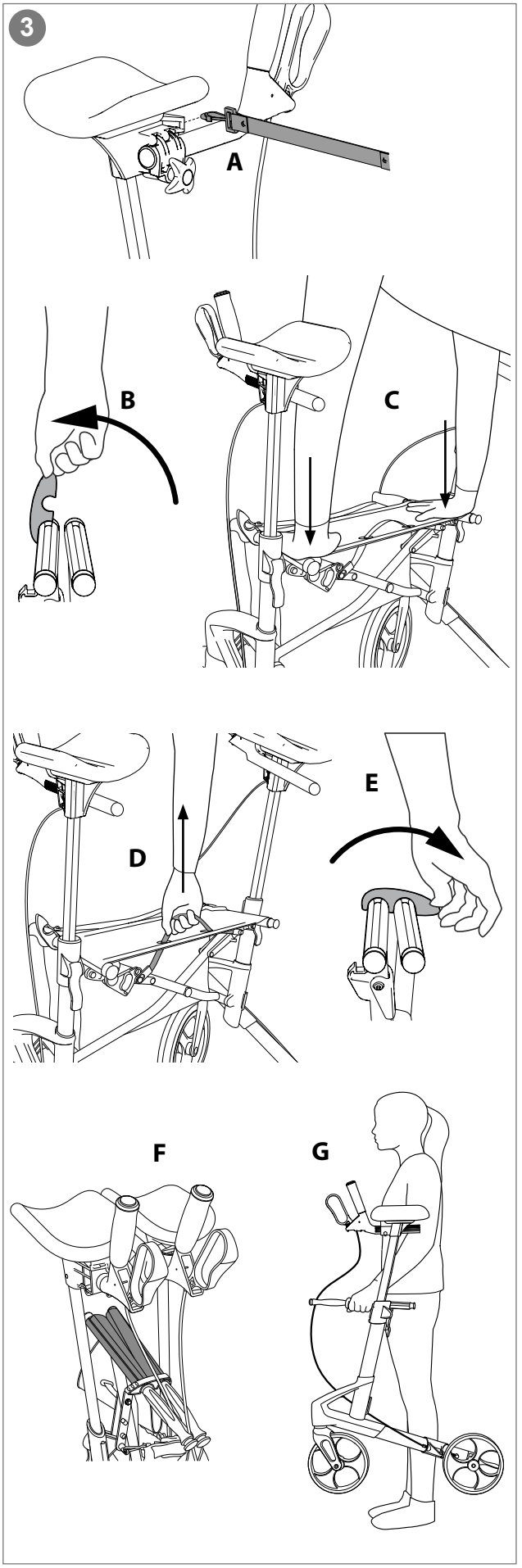
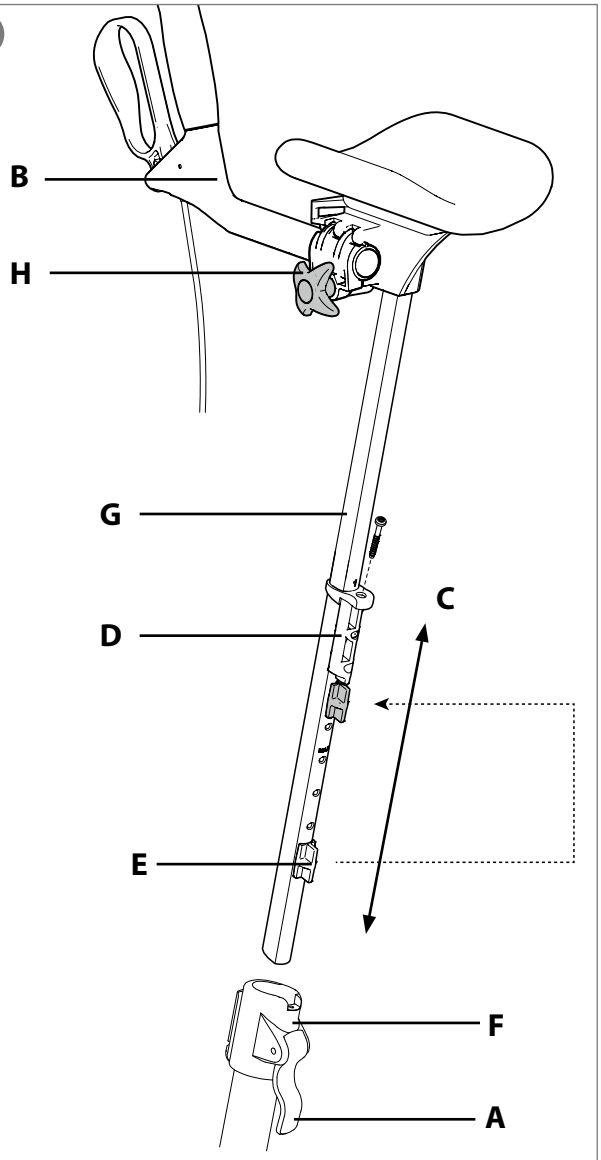
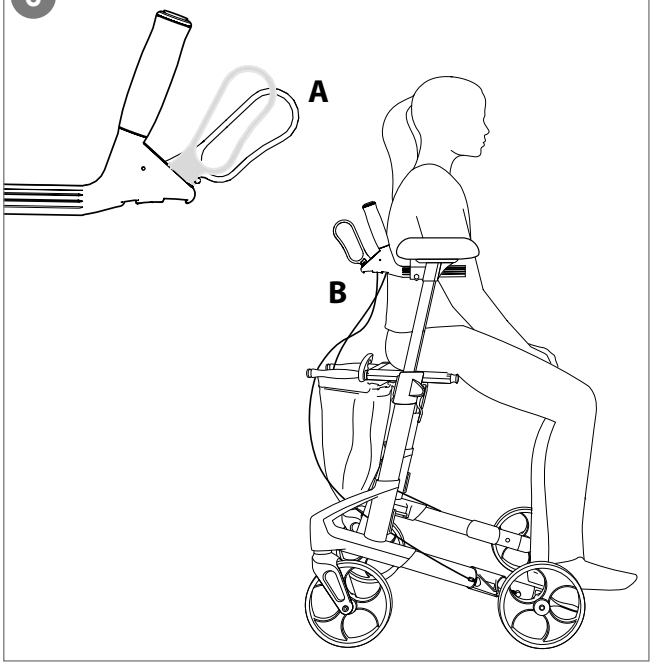
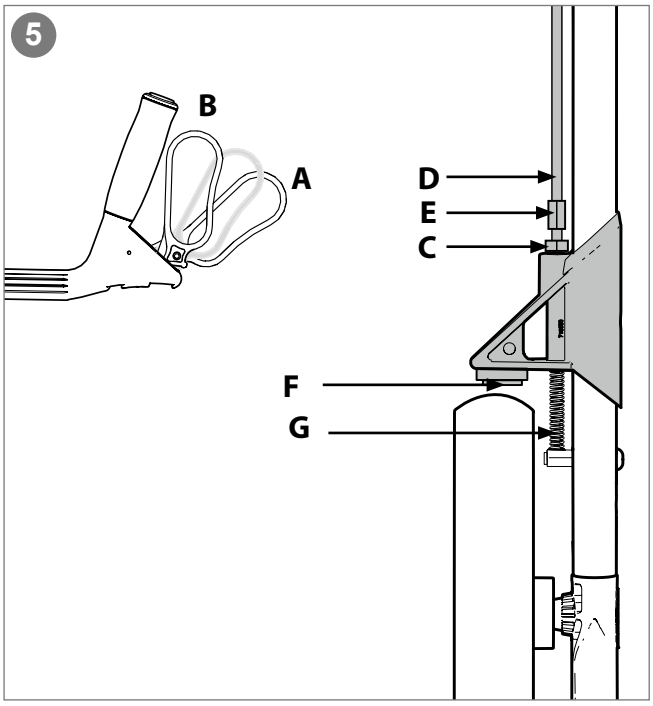
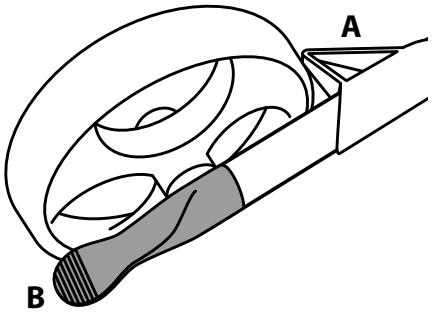
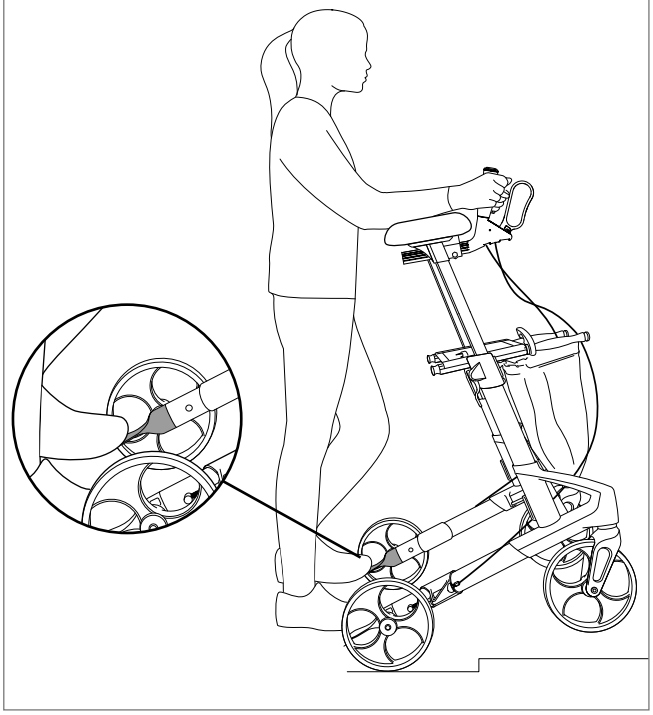
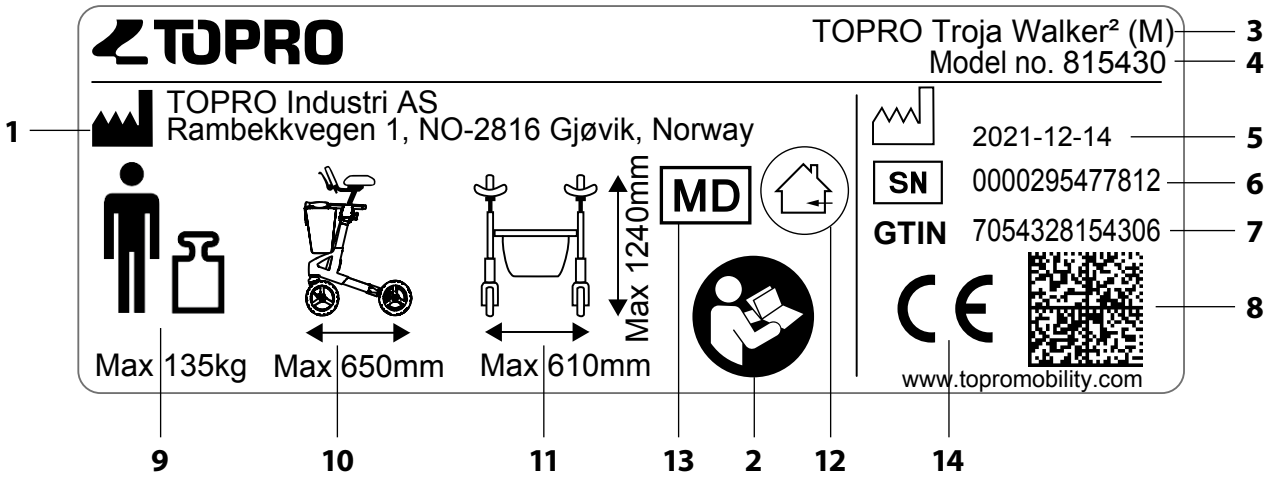
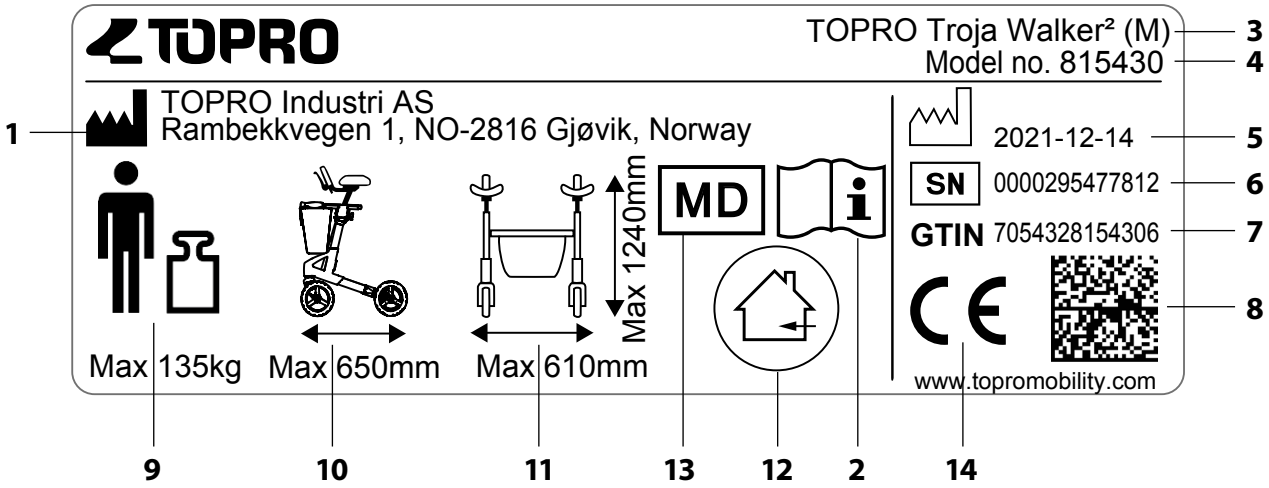
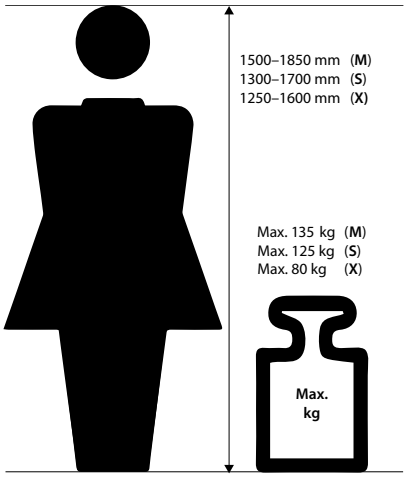
## Congratulations with your choice of a new TOPRO forearm walker

Every TOPRO forearm walker is developed, designed and manufactured in Norway. The product conforms to the Medical Device Regulation (EU) 2017/745. It is tested and approved according to EN ISO 11199-3:2005. On receipt of the goods, please check your product. If you have any queries, please contact your mobility shop or dealer immediately. If you have problems reading this user manual, an electronic version is available on **topromobility.com.**

## Intended use

The TOPRO Troja Walker² shall give support to users with reduced balance and/or reduced walking ability. It is designed to be pushed, not pulled. The integrated seat provides an opportunity for the user to rest if needed. The intended user is an adult. Height and weight restrictions apply. *Contraindications: The product is not suitable for persons with low strength in arms, with very poor balance, or with considerable cognitive disabilities. The walker is designed and approved for indoor use on plain foundation*.

To get to know the walker and learn how to use it properly, training by a professional is recommended.

## 0 The grey numbered circles in this user manual refers to the numbered information graphic and illustrations inside the cover

## 1 Get to know your TOPRO product

**A** Forearm support
**B** Brake lever
**C** Handle tube
**D** Brake cable
**E** Seat bars / carrying handle
**F** Seat with strap for folding mechanism
**G** Locking clamp
**H** Spak för höjdjustering
**I** Basket
**J** Side frame
**K** Brake block / Edge guard
**L** Tilting function
**M** Supporting cord
**N** Product label

## 2 Safety ⚠

When folding and unfolding the walker, be aware so you do not hurt your fingers.

Be aware not to put your fingers in open holes or between stationary parts.

Make sure that the levers for height adjustment are positioned downwards.

When stationary or parked, make sure the parking brakes on both sides are activated (A).

Do not use the walker on stairs. Be careful when passing thresholds etc.

Check the brakes before use every time, and that the walker locks in unfolded position.

Do not use the walker to transport heavy loads or persons.

Be aware that some parts of the walker might feel hot or cold when exposed to extreme temperatures.

Do not modify the product as this might put your safety at risk and the guarantee will become void.

Before sitting on the seat (C), make sure the parking brakes on both sides are activated (A).

To reduce the risk of falling, do not lean or bend sideways or over to the back of the walker whilst in the seated position (D). Do not tilt the walker too much sideways, to avoid the risk of falling. For the same reason, bags etc. must never be hung on the handles.

To give the right support and safety, walk straight and make sure the walker is kept as close to your body as possible (E). Walk between the rear wheels (F), not behind them (G).

Be aware that the walker may roll faster than you when walking downhill, which may put you in risk of falling. Activate the driving brakes to adjust the speed (B).

Be careful when wearing fluttering skirts or trousers as they could get caught in the wheels, which could lead to falling. Wear sturdy footwear to avoid injuries.

Do not use the walker on slippery surfaces.

Keep the walker away from fire and hot items.

Do not climb on the walker.

Art. No. User Manual: 104848 – Revision D, 2026-06

## 3 Unfolding / Folding

**Unfolding:** Loosen the strap between the forearm supports (A) (not applicable for the walker in size X). Release the lock clamp positioned at the seat bars (B). Press the seat bars down and outwards (C). The side frames will slide apart and you will hear a clicking sound. The walker is now locked in its open position. Fasten the strap (A).
**Folding:** Loosen the strap between the forearm supports (A). Pull the strap on the seat (D) and press the side frames together. Lock the clamp (E) over the seat bars.
**Carrying the walker:** The walker must not be carried by the seat strap or the brake cables. Fold the walker together as shown on the illustration (F). When folded, carry the walker by the seat bars (G).

## 4 Forearm supports / Handles / Memory function

**Height adjusting:** Pull the lever (A) upwards and adjust the forearm support tubes to a suitable height. Then move them slightly up and down (C) until you hear a click. Ensure that the lever is pointing downwards (A). Set both tubes to the same height.
**Recommended height:** Stand between the rear wheels. Relax your shoulders and hold your elbows in a 90-degree angle. The level of the forearm supports should be fixed at the level of your elbows. A professional should help you with the adjustment of the forearm supports and the handles, to help ensure correct posture.
**Adjusting the handles:** Pull the lever (A) upwards. You can now pull the handle (B) backwards or push it forwards. It is also possible to adjust the angle of the handle. Once the handle has been adjusted to the desired position, return the lever (A) to the downward position. Be careful when opening and closing the lever (A) to avoid injury to your fingers.
**Memory function:** Set the height of the forearm support tube. Count the visible holes on the rear of the forearm support tubes (G). Loosen the screw on the top of the attachment piece (D). Pull the lever (A) upwards and pull the forearm support tube (G) out of the frame. The attachment piece (D) will now follow the forearm support tube. Move the attachment piece (D) so that the number of holes above the attachment piece corresponds to the number of holes counted previously. Insert the red memory clip (E) into the first available hole below the attachment piece (D). Pull the lever (A) upwards and push the forearm support tube (G), together with the attachment piece (D), down into the side frame (F). Reattach the screw to the attachment piece (D). Pull the lever (A) upwards and pull the forearm support tube (G) upwards until it stops. Then release the lever. If the forearm support tubes are in the lowest position, they will automatically stop at the preset position when pulled upwards.

## 5 Brakes

Make sure the brakes are working before every trip. Each brake lever (A) works on each of the rear wheels. *Do not drive or push the walker while the parking brakes are activated.*
**Parking brakes:** Push both brake levers forward (A) until they lock. Be aware that the lever has some resistance before it stops and locks the rear wheel. Pull the brake levers backwards to release the parking brakes.
**Driving brakes:** Pull both brake levers towards you (B) to reduce speed. *Do not use continuously while the walker is in motion.*
**Adjusting / checking:** Loosen the nut (C). Tighten the brake cable (D) with the adjustment screw (E). Adjust both brake blocks (F) so the distance to the wheel surface is approximately 1 mm and then tighten the nut (C). Make sure the brakes are not adjusted too tight and pressing hard into the rubber of the wheels when the parking brake is activated.
**Maintaining:** Clean the brake block (F) and the brake spring (G) regularly. The wheels and brakes can be replaced if they are worn out or damaged.

## 6 Seat

The walker has a safe and convenient mesh seat, which provides view to the ground while walking. It can easily be wiped with a cloth if needed. The parking brakes must be activated before using the seat. Push both brake levers forward until they lock (A). Both rear wheels are now locked. Sit with your back towards the basket (B).

## 7 Edge guard / Tilt function

**Edge guard:** The brake blocks (A), which are in front of the rear wheels prevent the rear wheels from catching on corners, door frames etc.
**Tilt function:** To use the tilt function use the front of the foot and press down on the tilt pedal (B), which is beside the rear wheel. This helps to tilt the walker over thresholds, edges etc. Please note: Be careful when passing thresholds etc.

### Basket

The basket holds up to 5 kilos. It can be removed and washed gently at 60 °C.

### Transport

The walker should be transported in the folded and locked position (*see section Unfolding / Folding*). Ensure that the lock clamp is on, so that the walker cannot unfold unintentionally. Be careful when loading the walker into a vehicle, giving attention to securing it.

### Storage

The walker should be stored in an upright position. If stored outdoors, then under cover to protect it. Do not leave your walker unattended. Do not put heavy items on top of the walker during storage.

### Cleaning

The walker can be cleaned using a damp cloth. In case of mud and grime accumulation, a cloth with a mild washing detergent can be used followed by rinsing with lukewarm water. Do not use high pressure cleaners and avoid flushing directly against ball bearings (in wheels and front forks). Do not use abrasives or steam wash. Dry the walker using a soft cloth. When cleaning in automatic washing systems the following applies: *Washing; 60 °C, 10 min. (pressure maximum 4 bars, pH neutral) – drying; 70 °C, 10 min.*

## Disinfection

Disinfection must only be carried out by authorized personnel. Please contact your local dealer or mobility aid centre. During disinfection sufficient protective equipment must be worn, and the instructions of the disinfectant manufacturer must be followed. The surface of the walker should be cleaned with a soft cloth or sponge and an approved chemical disinfectant containing 70–80 % ethanol. We advise against disinfectant that contains chlorine or phenol. Do not use heat disinfection. *The manufacturer cannot be held responsible for any damage or injuries that might occur using harmful disinfectant/measures or disinfection carried out by unauthorized personnel.*

## Inspection / Maintenance / Re-use

It is recommended that inspection and maintenance are carried out regularly (frequency is determined by how and how often the walker is used). Check the following items: Frame, screws, tubes, handles, forearm supports, brakes, brake parts, wheels, seat and accessories. This also applies when the walker is to be prepared for re-use. Please see detailed instructions in this user manual about maintenance of brakes and cleaning/ disinfection. The recommended maintenance is not a requirement, and no preventive maintenance is required provided that the walker is used as intended according to this user manual.

## Material / Recycling

The walker is made of plastic-coated aluminium tubes, plastic couplings and material made of plastic and polyester. There are ball bearings in all wheels and forks. Most parts can be recycled. Dispose the walker and its packaging according to the applicable regulations in your country. For information, please contact your public authorities.

## Warranty

TOPRO Troja Walker² is guaranteed free of faults and defects for 7 years. Damage caused by incorrect use or parts that are exposed to natural wear and tear (e.g. brake block, brake cable, wheel, basket, seat and grip) are exempt from the 7-year warranty. For repairs during the warranty period, please contact your local mobility shop or dealer. The warranty is void if unauthorised spare parts or accessories have been or are being used on the product. The expected life span of the product is estimated to be 10 years if it is used in the correct manner according to this user manual and its safety and maintenance instructions.

## 8 Product label

1 - Manufacturer
2 - Consult Instructions for use (blue icon on the new label)
3 - Model name
4 - Model number
5 - Date of manufacture
6 - Serial number
7 - Global Trade Item Number
8 - GS1 DataMatrix
9 - Maximum user weight
10 - Maximum length of walker
11 - Maximum width / height of walker
12 - Intended for indoor use
13 - Medical Device
14 - CE-marking (This product conforms to Medical Device Regulation (EU) 2017/745)

## Accessories

| | TOPRO Item No. |
|---|---|
| Back support | **814765** |
| Tray (NB: not for size X) | **814728** |
| Anti-slip mat for tray (NB: not for size X) | **814059** |
| Crutch holder | **814040** |
| Bottle holder | **814043** |
| Stabilizer bar (size M / S) | **814729** |
| Stabilizer bar (size X) | **814730** |
| One-handed brake | **814026** |
| Guiding handle for the visually impaired | **814014** |
| Rear bag with zipper (NB: not for size X) | **814046** |
| Name tag | **814024** |
| Roll Resistance System (drag brake), pair | **814032** |
| Basket for oxygen bottle (only compressed, not liquid) | **814009** |

Accessories mounted on the walker may influence stability, we therefore advise you to use them with care. Accessories can be ordered separately to enable a TOPRO Troja Walker² to be customised to individual needs. Contact your mobility shop, dealer or TOPRO for an updated overview of available accessories or visit our homepage **topromobility.com.**

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